
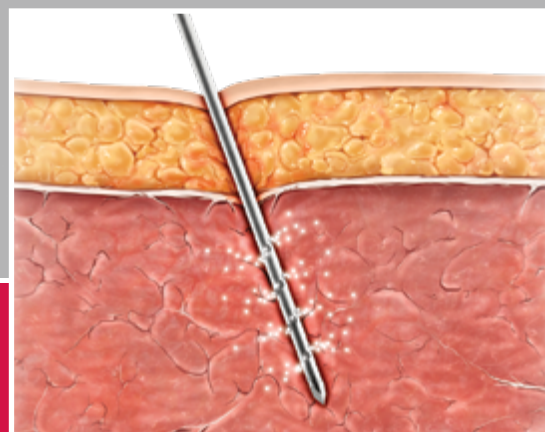


# The Difference in Therapeutics Delivery

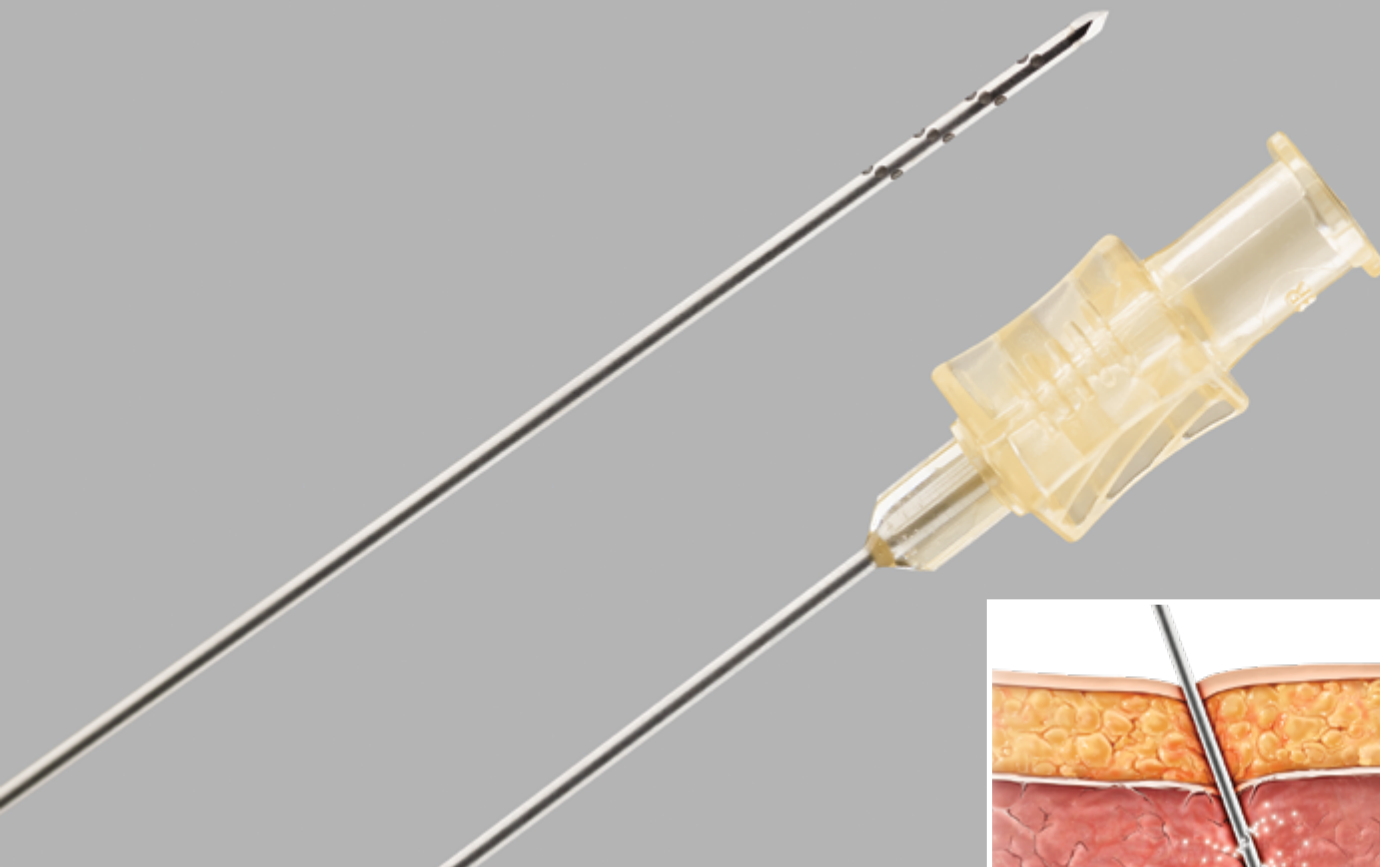

## ProFusion™ THERAPEUTIC INFUSION NEEDLE

The ProFusion Therapeutic Infusion Needle is used for intramuscular or subdermal infusion of diagnostic and therapeutic agents. The difference is evident in:

- The 1cm distal infusion segment has sideports arranged in a spiral to help uniformly deliver diagnostic and therapeutic agents
- The stainless steel needle has a trocar tip to aid with precise placement
- The needle sideports are highly visible under ultrasound to allow visualization for targeted agent delivery

For more information about the difference in intramuscular or subdermal infusion of diagnostic and therapeutic agents, ask your Cook Regentec representative.




# ProFusion™ Therapeutic Infusion Needle

Order Number	Reference Part Number	Needle Gage	Needle Length (cm)	Sideport Diameter (mm)
G47570	PTIN-19-5.0	19	5.0	0.56
G47571	PTIN-19-7.5	19	7.5	0.56
G47572	PTIN-19-10.0	19	10.0	0.56
*G47573	PTIN-19-15.0	19	15.0	0.56
G47574	PTIN-21-5.0	21	5.0	0.46
G47575	PTIN-21-7.5	21	7.5	0.46
G47576	PTIN-21-10.0	21	10.0	0.46
*G47577	PTIN-21-15.0	21	15.0	0.46
G47578	PTIN-25-2.5	25	2.5	0.18
G47579	PTIN-25-5.0	25	5.0	0.18
*G47580	PTIN-25-7.5	25	7.5	0.18

\*Coming soon.

The ProFusion Therapeutic Infusion Needle, a Cook Regentec tool to help build your therapies.

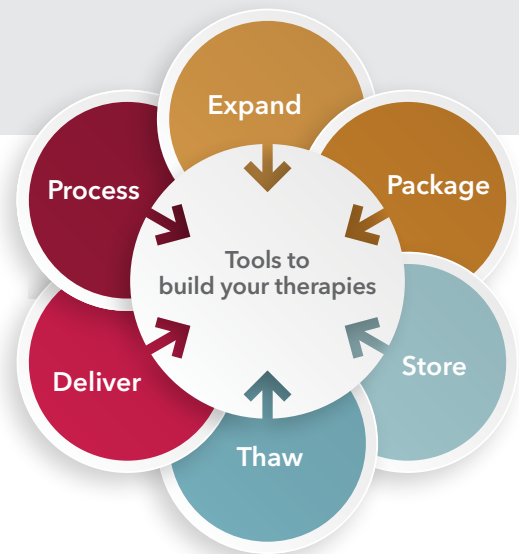

From concept to commercialization, Cook Regentec creates tools to help you build your therapies.

### Customer Service

Phone: +1 317.917.3450  
 Toll-free in the United States: +1 800.265.0945  
 E-mail: [customer-service@cookregentec.com](mailto:customer-service@cookregentec.com)

Cook Regentec  
 1102 Indiana Avenue  
 Indianapolis, IN 46202 U.S.A.

This product may not be available in all regulatory jurisdictions.



[www.cookregentec.com](http://www.cookregentec.com)